


AID2GO: User Manual/Quick Start Guide

Features:

- 16x2 LCD Display
 - Speaker
 - 4 Buttons for User Input
 - Adjustable Time Settings
 - Store up to 99 User Alerts
 - Enable/Disable Alerts at any time
 - 9V Power Supply
-

Button Reference:

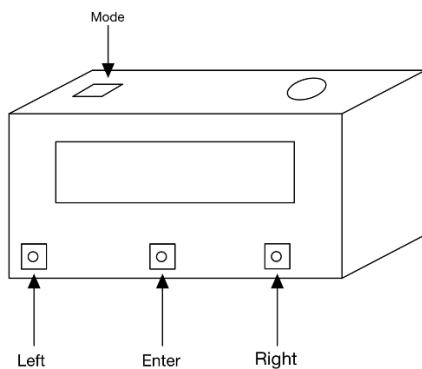

Getting Started:

Once you receive your AID2GO, you first want to **set the time** and **some alerts**.

Your device will already be in the **Main Menu**



Press **enter**

Your device will now be in the **time set** mode



From here, you can use the **enter** button to change the value you are editing (year, month, etc.)



Use the **left/right** buttons to adjust this value. Once you are done setting the time, press the **mode** button to go to the **Alerts Menu**.



Use the **left/right** buttons to change the reoccurring event that you are editing.



When you are on the event you want to edit, press **enter** to begin changing the properties.



Make sure to set the **EN (enable)** value to **ON** if you want the alert to trigger.



When you are done, press the **Mode** button to return to the **Main Menu**.

Receiving Alerts:

When an alert is triggered, a screen displaying the amount of missed alerts, the type of alert, the alert number, and the time it was missed will be displayed. A tune will also play during this screen.



To “acknowledge” this alert, press the **enter** button, and the AID2GO will return to the main menu and stop playing the tune. You do **not** have to re-enable an alert after you clear it, it will repeat every week unless manually disabled in the menu.

Time Settings:

- Year (YYYY)
- Month (MM)
- Day (dd)
- Hour (hh, 24h set, 12h display)
- Minute (mm)

Reoccurring Alert Properties:

- Alert Number (0-99)
- Day of Week (Sun-Sat)
- Hour (hh, 24h)
- Minute(mm)
- Type (Alert, Doctor, Medicine)
- Enable (On/Off)

Having Troubles?

Call our support desk at (123) 456-7890 for advanced troubleshooting and general assistance.